


HyBase 1000


Mechanical & Hydraulic Operating Table

Technical Specifications


Dimensions	
Full length	2000mm
Height without padding	780mm-1020mm
Width without side rail	520mm
Width with side rail	590mm

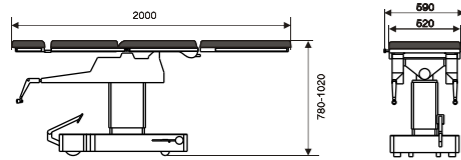


Trendelenburg position	
Trendelenburg	Max. 22°
Reverse trendelenburg	Max. 22°

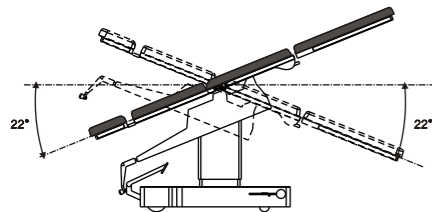


Leg plates	
Leg plates up	Max. 15°
Leg plates down	Max. 85°

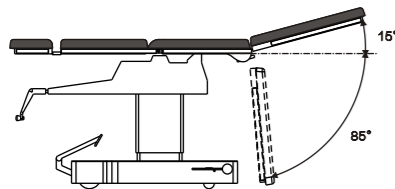


Reflex Position	
Reflex Position	105°

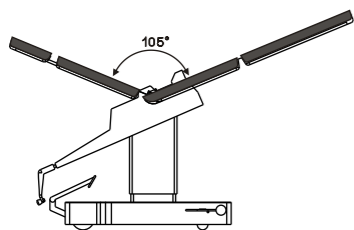


Lateral tilt	
Left tilt	Max. 18°
Right tilt	Max. 18°

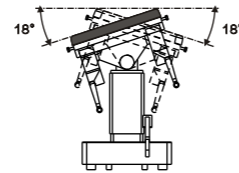


Back plate	
Back plate up	Max. 75°
Back plate down	Max. 15°

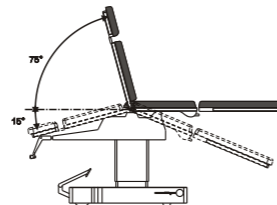

Head plate	
Head plate normal travel	Max. 45°
Head plate extended travel	Max. 70°

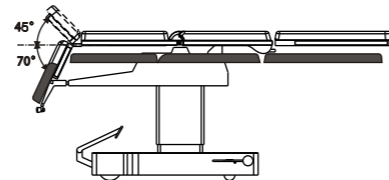

Flex Position	
Flex Position	225°

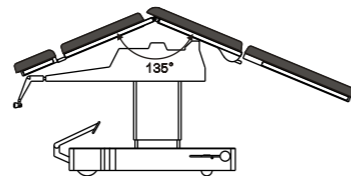

HyBase 1000 Operating Table

NEW



DISTRIBUTOR:



Mindray is the registered trademark owned by Shenzhen Mindray Bio-medical Electronics Co., LTD. Hybase is the registered trademark owned by Nanjing Mindray Bio-medical Electronics Co., LTD © 2009-2010 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights Reserved. P/N:ENG-HyBase 1000-420285-20100525

mindray

Mindray is listed on the NYSE under the symbol "MR"
Mindray Building, Keji 12th Road South, High-tech Industrial Park,
Nanshan, Shenzhen 518057, P.R. China
Tel: +86 755 26582888 Fax: +86 755 26582680
E-mail: intl-market@mindray.com Website: www.mindray.com

mindray
healthcare within reach

Mindray OR Total Solution

All the advanced and reliable facilities from Mindray technology, are perfectly matched to give you total support!



HyBase 1000

Mechanical & Hydraulic Operating Table



- Persistence in Mindray classic quality and expertise apperception in operating table, perfect flexibility and agility, make it totally cater to different kinds of operating environment requirements, and make it easy for your operation.
- Ingenious gear transmission design and excellent gear match make all movements to be done in a blink without any trouble.
- With sufficient ergonomic consideration, the base is T-shaped and capable of providing maximum space for the surgeons' feet so that they can stand through long time operations comfortably.



Table up and down foot pedal



Table locking & release foot pendal



HyBase 1000

Mechanical & Hydraulic Operating Table



Table caster



Table break system



Flex position



Table up



Table down



Re-flex position




Trendelenburg position



Reverse Trendelenburg position



Leg plate up and down, split 180°

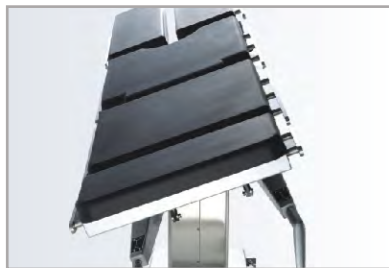

Turn left



Turn left



Stainless steel rail & antistatic pads



Back plate up



Back plate down



Head rest up



Head rest down



180mm-high body elevator for kidney and gallbladder operations



HyBase 1000

Mechanical & Hydraulic Operating Table



High level stainless steel structure for 135kgs weight capacity

- The main body of the operating table is made of stainless steel, which offers excellent corrosion resistance, thus making it easy to clean.
- Designed conforming to DIN standards, the rails are compatible with all series accessories application.



ENT & Neurosurgery position



Beach chair position

- The removable pad is made of PU, a special antistatic and antimicrobial material, and is easy to clean and comfortable to the patients.
- The robust stainless steel frame provides absolute stability and a safe environment for the patients.




High quality builds up high confidence

HyBase 1000

Mechanical & Hydraulic Operating Table



Normal position for Spinal operation

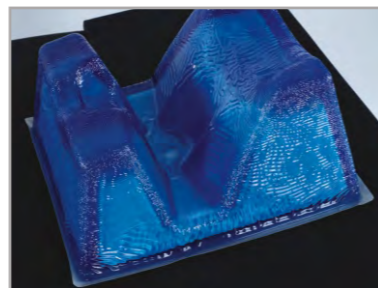

115-005830-00
Prone head support



115-005885-00 Big horseshoe pillow
115-005887-00 Middle horseshoe pillow
115-005889-00 Small horseshoe pillow



Reverse position for Spinal operation



115-005788-00
Intervertebral disc operation pad



Lithotomy position for Gynecology & Urology



115-005910-00
Pad for lateral position



115-005814-00
Heel support pad



Lateral position



Lateral position